

ThermalSeal® RTS™


Features a pressure-activated, non-tacky, silicone adhesive

i Application:

qPCR, long-term storage, protein crystallography & HPLC

🌡️ Temperature Range:

-80°C to +105°C

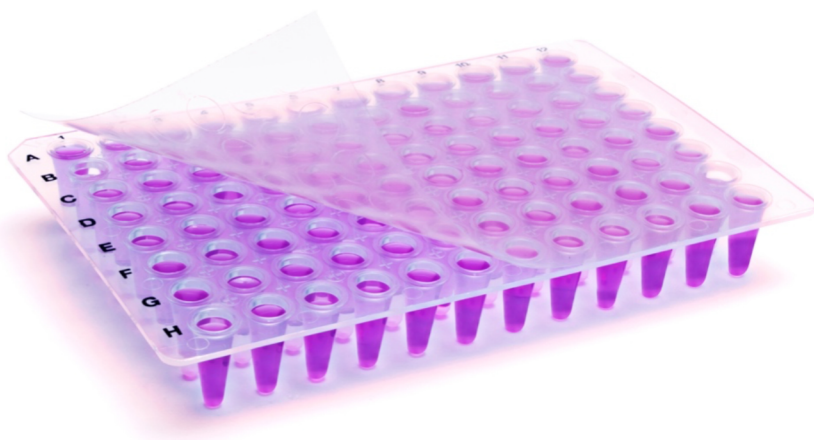

CATALOG #	DESCRIPTION	QUANTITY
TSS-RTQ-100	ThermalSeal RTS	100 Sheets
TSS-RTQ-25	ThermalSeal RTS	25 Sheets
TSS-RTQS-50	ThermalSeal RTS, Sterile	50 Sheets

Features and Benefits

- Certified DNase- & RNase- free
- Handling and plate positioning is simplified as a result of the non-tacky adhesive
- Chemically resistant to DMSO, acetonitrile, and methanol
- Ideal for use with flat-top, raised-rim and non-skirted PCR plates
- High optical clarity and low auto-fluorescence
- The silicone adhesive is non-tacky until pressure is applied to the sealing surface which then creates a secure bond around the wells and perimeter of the plate
- Ideal for dry samples (e.g. soil, seeds, insects)
- A roller or sealing paddle is highly recommended for use with this product to ensure consistent and

Bundling Opportunities

- PCR plates
- Protein crystallography plates
- Deep well plates
- PCR reagents and buffers
- Thermal cyclers
- Pipettors and pipet tips
- Reservoirs

